


Installation instructions



Hot Tracker



For your surgical trays.

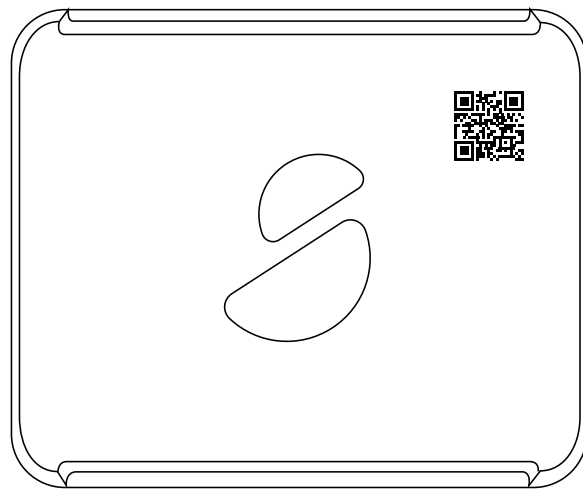


Handle the tracker with care when
attaching it to the tray to avoid injury



Supplied Tools

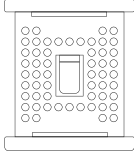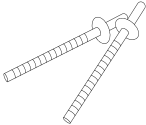	Hot Tracker 78*65*33.5 mm 3.07*2.56*1.32 inch
	Bracket The size of the bracket is adapted to the tracker's dimension.
	x2 rivets Ø 4 x L6 to be tightened with a standard rivet clamp. Rivets are supplied by SMADE

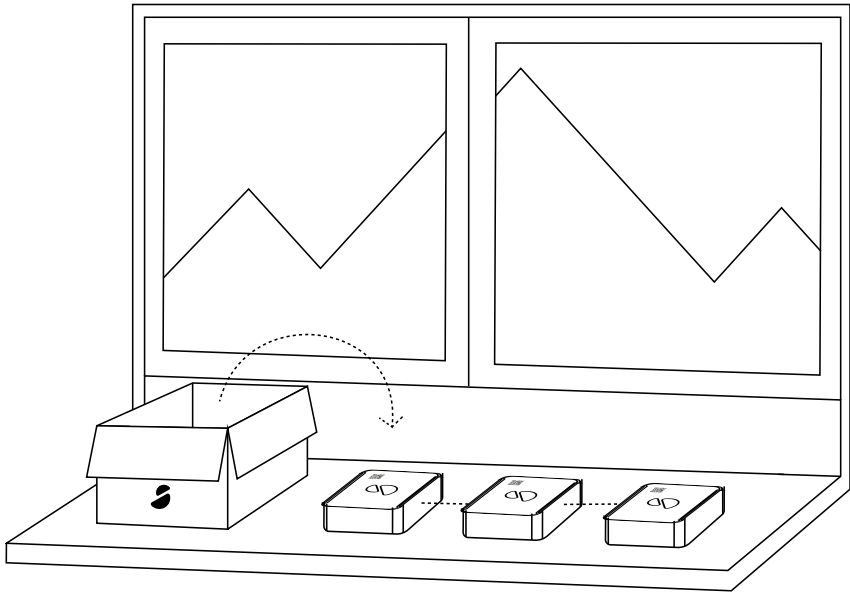


For optimal performance, remove the trackers from their packaging, place them flat (avoid stacking), and keep them near windows or glass doors — not in basements or enclosed areas

1

Identify a free space in the tray to install the tracker.


The tracker should ideally be installed on the side walls but can alternatively be located on the bottom wall.

Make sure not to place the tracker in the space where an instrument may be inserted.

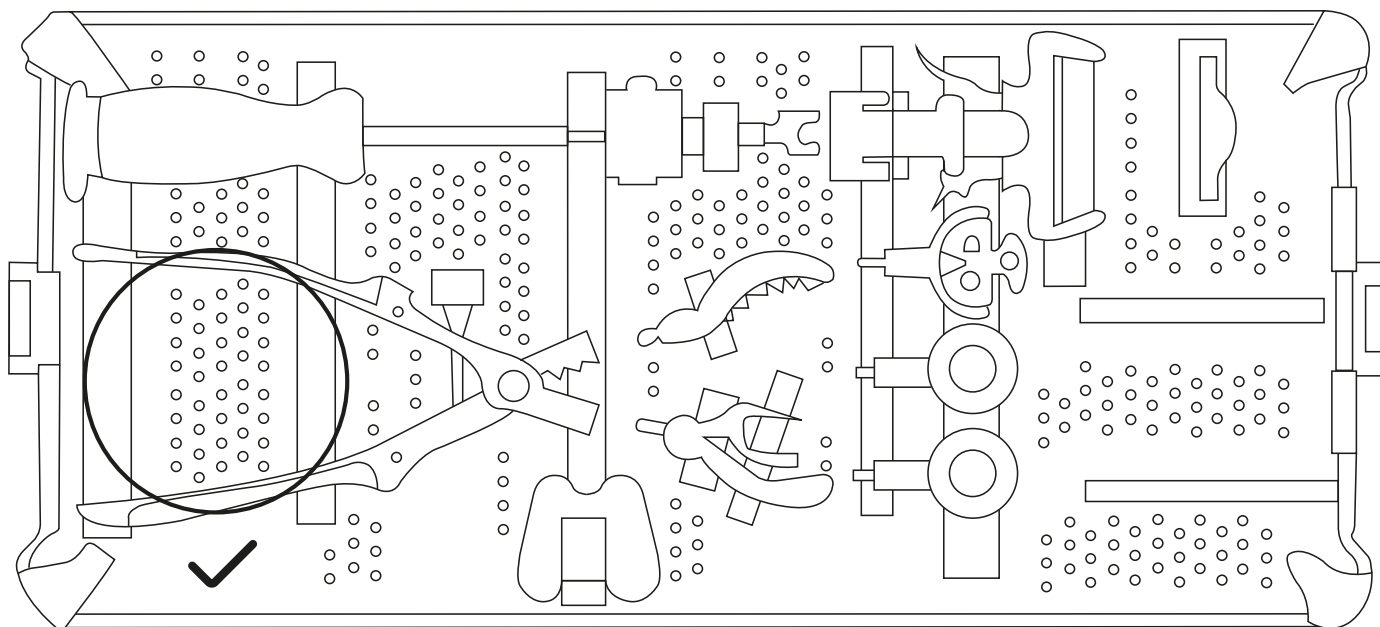


2


Once location identified, place the bracket to check perforation alignment between the bracket and the tray.

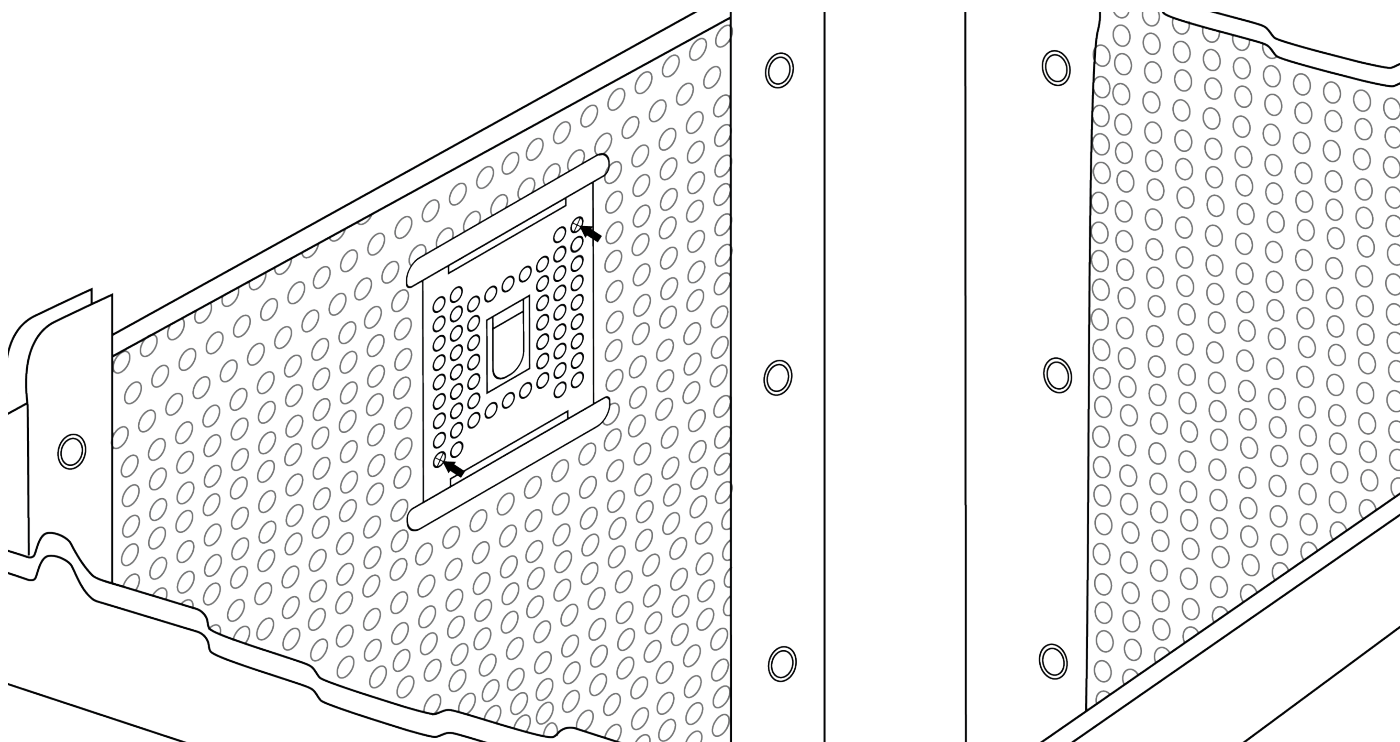

- 3** Place 2 rivets to secure the bracket to the tray wall.

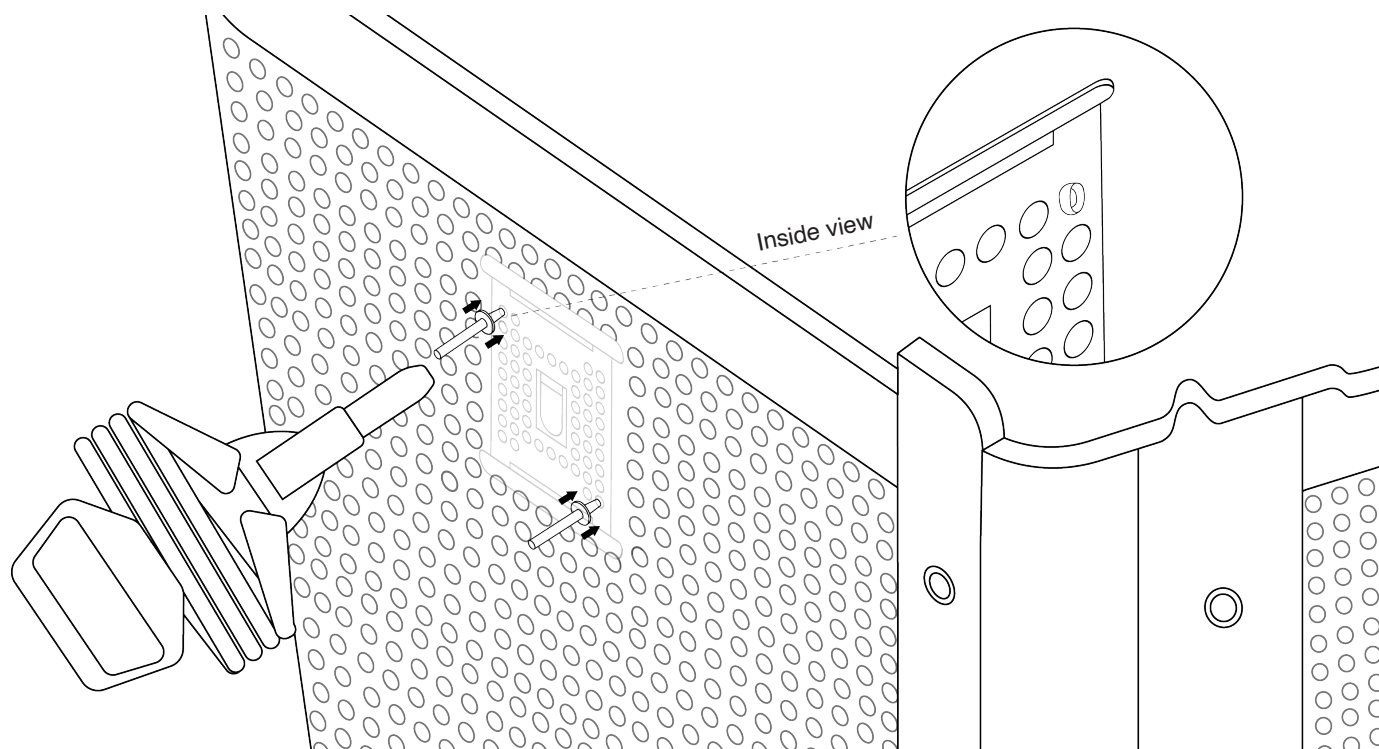

- 4** Slide the tracker's bottom end into the bracket. Press down on the top end of the tracker to securely snap into place inside the bracket.

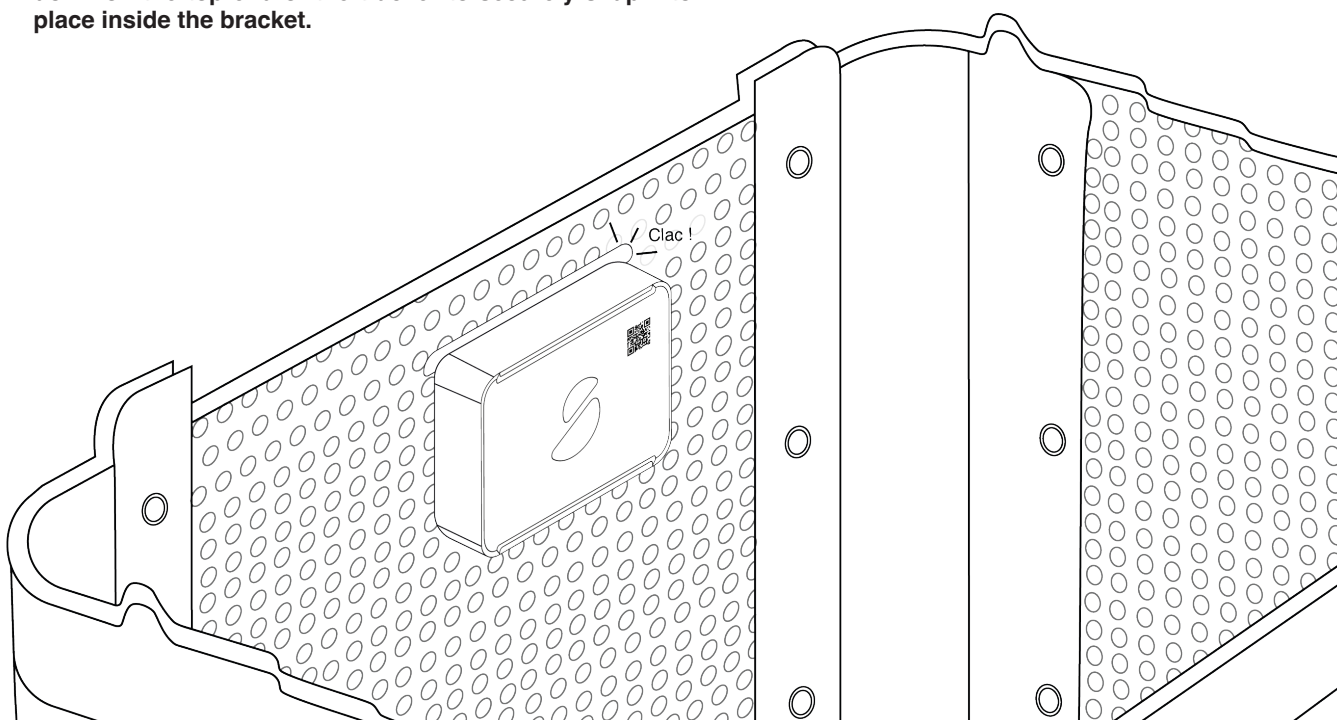


- ✓ To remove the tracker, open the bracket slightly (you may need to use a small tool to do so).
- ✓ Tracker can be removed and repositioned if needed and when required.
- ✓ The bracket may need to be adapted to maintain appropriate pressure on the tracker.

 **Tracker Attachment — Watch and Learn!**



Hot Tracker

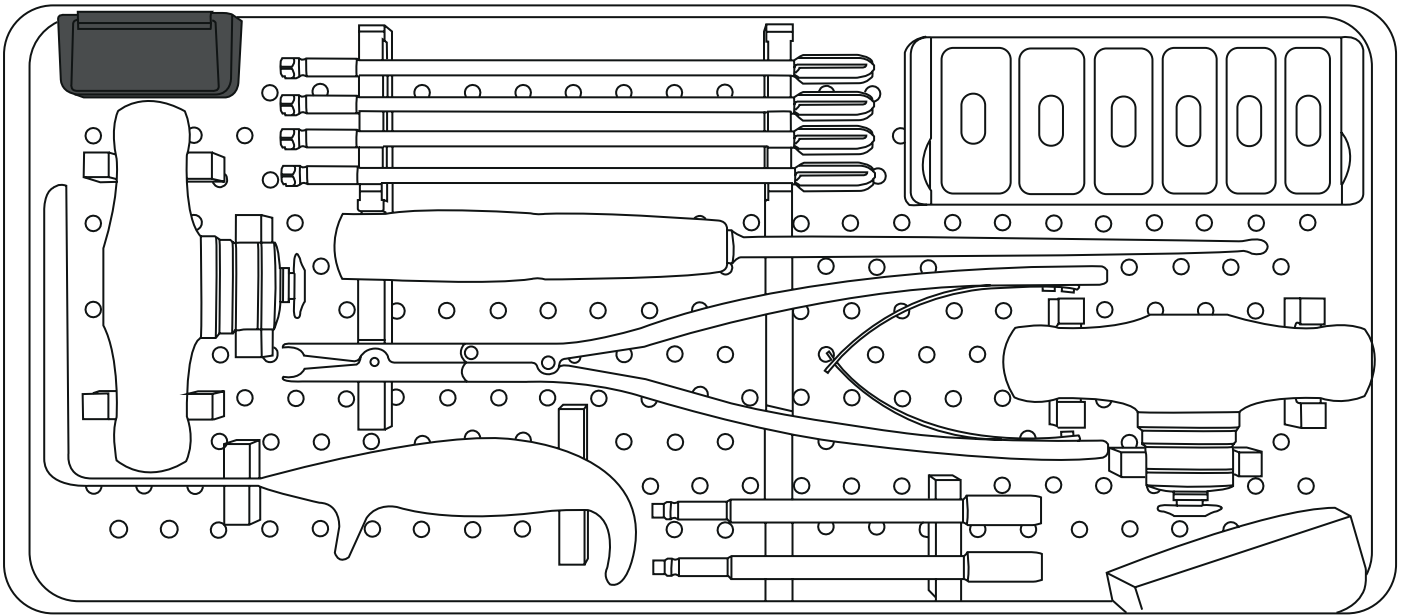

If you have any specific installation-related needs, rest assured that we're more than happy to assist you. Simply share photos or videos with us at mike@smade.io or we can schedule a personalized meeting to help you find the best solution.